




JDENTALCARE
just smile

JD ICON PLUS T

The power of T design

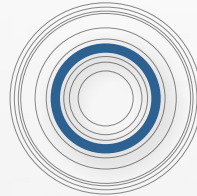


conical connection



Made in Italy

PRODUCT FEATURES

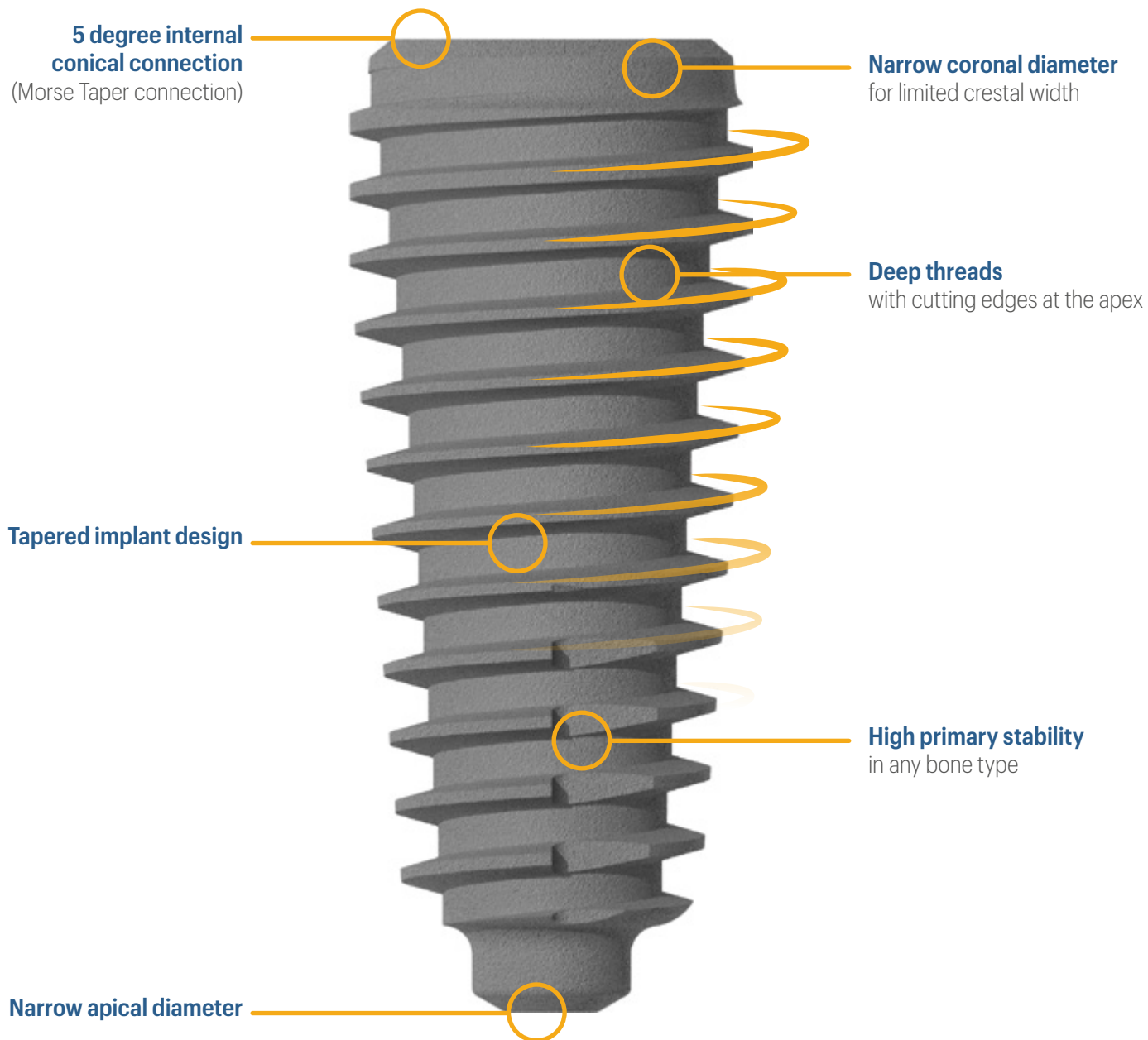


PRODUCT SPECIFICATIONS

The JDIcon Plus T dental implants are available in the diameters and lengths shown in the following table:

IMPLANT DIAMETER	TIP DIAMETER	ABUTMENT INTERFACE	PLATFORM DIAMETER	BEVEL DIAMETER	WIDEST THREAD DIAMETER	LENGTH			
						8,5	10	11,5	
Ø 3,5	1,6	3,1	3,4	3,8	4,0		10	11,5	13
Ø 4,0	1,75	3,1	3,6	4,0	4,4	8,5	10	11,5	13
Ø 4,5	1,75	3,1	3,6	4,0	4,9	8,5	10	11,5	13
Ø 5,0	2,2	3,1	4,0	4,5	5,5	8,5	10	11,5	13
Ø 6,0	3,0	3,1	5,0	5,5	6,4	8,5	10	11,5	

Note: all measurements in mm

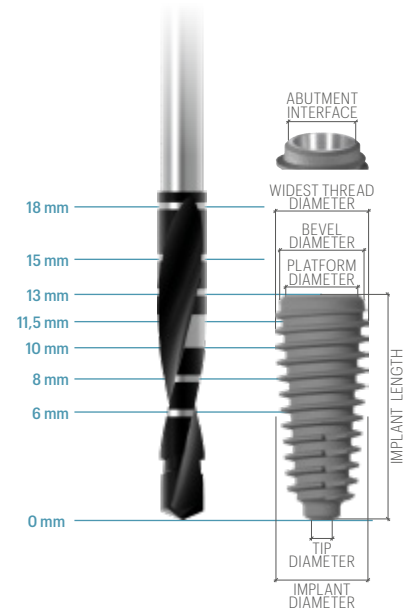


PRODUCT CATALOGUE

Note: all measurements of diameters and lengths are in mm

Implants:

IM37100-0	JDIcon Plus Ø 3.5 L 10
IM37115-0	JDIcon Plus Ø 3.5 L 11.5
IM37130-0	JDIcon Plus Ø 3.5 L 13
ICT40080.	JDIcon Plus Ø 4.0 L 8.5
ICT40100.	JDIcon Plus Ø 4.0 L 10
ICT40115.	JDIcon Plus Ø 4.0 L 11.5
ICT40130.	JDIcon Plus Ø 4.0 L 13
ICT45080.	JDIcon Plus Ø 4.5 L 8.5
ICT45100.	JDIcon Plus Ø 4.5 L 10
ICT45115.	JDIcon Plus Ø 4.5 L 11.5
ICT45130.	JDIcon Plus Ø 4.5 L 13
ICT50080.	JDIcon Plus Ø 5.0 L 8.5
ICT50100.	JDIcon Plus Ø 5.0 L 10
ICT50115.	JDIcon Plus Ø 5.0 L 11.5
ICT50130.	JDIcon Plus Ø 5.0 L 13
ICT55080.	JDIcon Plus Ø 6.0 L 8.5
ICT55100.	JDIcon Plus Ø 6.0 L 10.0
ICT55115.	JDIcon Plus Ø 6.0 L 11.5

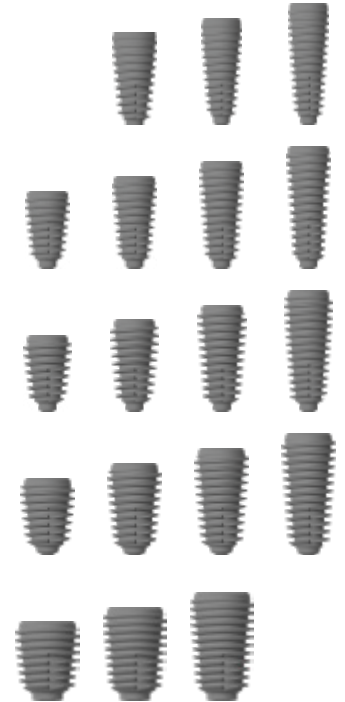

Note: All the prosthetic components are the same as JDIcon Plus implant line

DRILLING PROTOCOLS

It is recommended to adhere to the indications of the following drilling sequence to ensure optimal primary stability of the implants.

IMPLANT DIAMETER	HEALED BONE		POST EXTRACTIVE BONE	
	SOFT BONE	MEDIUM-DENSE BONE	SOFT BONE	MEDIUM-DENSE BONE
site preparation in maxilla				
Ø 3,5	2,0	2,0	2,0	2,0
	2,4 2,8 up to the 1 st laser mark L 6mm	2,4 2,8 3,2 up to the 1 st laser mark L 6mm	2,4 2,8 3,2 up to the 1 st laser mark L 6mm	2,4 2,8 2,8 at the entrance
Ø 4,0	2,0	2,0	2,0	2,0
	2,4 2,8 3,2 at the entrance	2,4 2,8 3,2 up to the 1 st laser mark L 6-8mm 3,6 at the entrance	2,4 2,8 3,2 at the entrance	2,4 2,8 3,2 at the entrance
Ø 4,5	2,0	2,0	2,0	2,0
	2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm	2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm 4,0 at the entrance	2,4 2,8 3,2 3,6 at the entrance up to the 1 st laser mark L 6mm	2,4 2,8 3,2 3,6 at the entrance up to the 1 st laser mark L 6mm
Ø 5,0	2,0	2,0	2,0	2,0
	2,4 2,8 3,2 3,6 4,0 at the entrance	2,4 2,8 3,2 3,6 4,0 up to the 1 st laser mark L 6mm 4,8 at the entrance	2,4 2,8 3,2 3,6 at the entrance	2,4 2,8 3,2 3,6 up to the 1 st laser mark L 6mm 4,0 at the entrance
Ø 6,0	2,0	2,0	2,0	2,0
	2,4 2,8 3,2 3,6 4,0 4,8 up to the 1 st laser mark L 6mm	2,4 2,8 3,2 3,6 4,0 4,8 up to the 1 st laser mark L 6mm	2,4 2,8 3,2 3,6 4,8 at the entrance	2,4 2,8 3,2 3,6 4,0 at the entrance 4,8 up to the 1 st laser mark L 6mm
site preparation in mandible				
Ø 3,5	2,0	2,0	2,0	2,0
	2,4 2,8 3,2 up to the 1 st laser mark L 6mm	2,4 2,8 3,2 up to the 2 nd laser mark L 8mm 3,6 up to the 1 st laser mark L 6mm	2,4 2,8	2,4 2,8 3,2 at the entrance
Ø 4,0	2,0	2,0	2,0	2,0
	2,4 2,8 3,2 3,6 up to the 2 nd laser mark L 8mm	2,4 2,8 3,2 3,6 4,0 up to the 2 nd laser mark L 8mm 4,4 up to the 1 st laser mark L 6mm	2,4 2,8 3,2 at the entrance	2,4 2,8 3,2 3,6 at the entrance
Ø 4,5	2,0	2,0	2,0	2,0
	2,4 2,8 3,2 3,6 at the entrance 4,0 at the entrance	2,4 2,8 3,2 3,6 4,0 4,4 4,8 at the entrance	2,4 2,8 3,2 3,6 at the entrance	2,4 2,8 3,2 3,6 4,0 at the entrance
Ø 5,0	2,0	2,0	2,0	2,0
	2,4 2,8 3,2 3,6 at the entrance 4,0 at the entrance	2,4 2,8 3,2 3,6 4,0 4,4 4,8 at the entrance	2,4 2,8 3,2 3,6 at the entrance	2,4 2,8 3,2 3,6 4,0 at the entrance
Ø 6,0	2,0	2,0	2,0	2,0
	2,4 2,8 3,2 3,6 4,0 4,8 at the entrance	2,4 2,8 3,2 3,6 4,0 4,4 4,8	2,4 2,8 3,2 3,6 4,8 at the entrance	2,4 2,8 3,2 3,6 4,0 4,8 up to the 1 st laser mark L 6mm



Note: all measurements in mm

JDentalCare, the JDentalCare logo and other trademarks mentioned in this document are, excluding specifically mentioned exceptions, trademarks property of JDentalCare.

The images of the products represented in this brochure are not necessarily to scale.

Note: Not all products manufactured by JDentalCare and represented in this brochure have been regulatory cleared/released for sale in all countries and/or have been licensed in accordance with FDA Rules and Regulations.

For more information about current product assortment and availability, please contact your local JDentalCare distributor. For prescription only.

Caution: Federal (US) law restricts this device to sale by or on the other for a licensed dentist.

See instructions for use for more information such as indications, warnings, precaution and contraindications.